


## THE AUTOINJECTOR SYSTEM

This guide is not intended to replace or supersede the advice of your Doctor. Always consult your doctor if you have questions or concerns involving penile injection therapy.

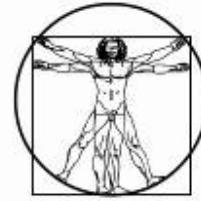

**MENS CLINIC INTERNATIONAL**  
SPECIALISING IN MALE SEXUAL HEALTH

## ABOUT YOUR NEW DEVICE


Penile injection therapy is a proven effective method of treatment for Erectile dysfunction.

The Autoinjector System is specially designed to help with the self-injection of medication into the corpora cavernosa of the penis.

This exclusive design of injection tip helps reduce sensitivity at the injection site. And the 'jet action' release provides a quick, smooth injection. We think you will be surprised and pleased with the speed, accuracy and overall comfort of penile injections with the Autoinjector System.

## AUTOINJECTOR COMPONENTS

Your Autoinjector kit contains:  
one injection aid device with auto-release button  
and one plunger cap.

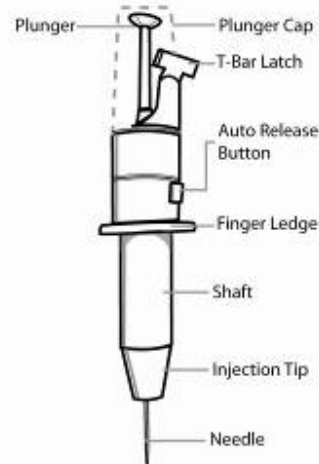

## USING YOUR AUTOINJECTOR

Before using your Autoinjector for the first time, you may wish to practice the application several times without a syringe, to become familiar and more comfortable with the technique.

Always wash your hands thoroughly before starting the self-injection procedure. Next check the medication vial for the expiration date, and discard if expired. Wipe the medication vial with an alcohol swab before proceeding.

## FILLING THE SYRINGE

Remove the needle cap from the syringe. Fill your syringe with a volume of air equal to the amount of medication that you will inject.

Pierce the rubber stopper of the medication vial and inject the air from the syringe into the bottle

**(See Figure 1).** *IMPORTANT: BE SURE THAT NO BUBBLES ARE IN THE FLUID PRIOR TO INJECTION.*

Then withdraw the amount of fluid prescribed by your doctor. Replace the needle cap over the syringe and snap into place.

## LOADING THE AUTOINJECTOR

Depress the T-bar to lock the syringe firmly in place (See Figure 2).

Remove the plunger cap, pull upward and back on the T-bar latch to open and cock the device. Put the filled syringe with the needle cap into the opening (See Figure 3).

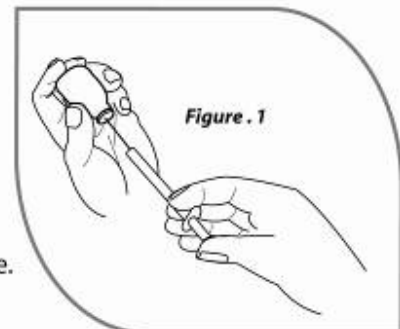


Figure . 1

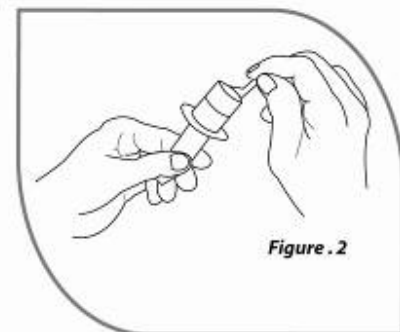


Figure . 2

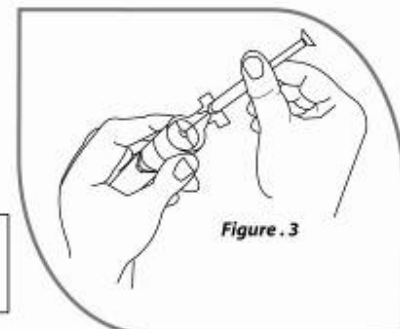


Figure . 3



**PROLONGED ERECTION!**  
SHOULD YOUR ERECTION LAST LONGER THAN FOUR (4) HOURS YOUR DOCTOR MUST BE CONTACTED IMMEDIATELY ON THE EMERGENCY NUMBER PROVIDED OR CONTACT 083 651 8672 (EMERGENCY ONLY AND IF DOCTOR IS UNAVAILABLE).

## THE SELF-INJECTION PROCEDURE

After filling the syringe and loading your Autoinjector, remove the needle cap. You are now ready to inject your prescribed dosage into the penis to create an erection.

Injections into the penis can be done from the right or left side, but always inject at either the **1:00** or **11:00** o'clock position. Avoid the top midline area, since this area contains several large veins. Even though you inject on just one side, medication still reaches the other side, since there are blood vessels connecting the two corpora of the penis. Avoid the urethra on the underside of the penis. (See Figure 4).

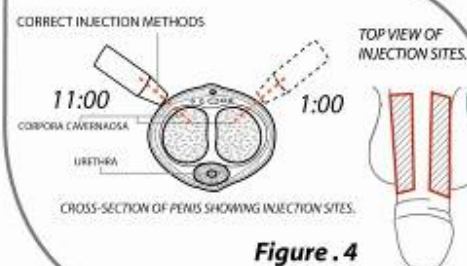


Figure . 4

### Step 1:

Wipe the injection site with an alcohol swab. Hold the Autoinjector by placing your thumb on the auto-release button and resting your middle finger on the finger ledge opposite your thumb (See Figure 5). This leaves your index finger free to push in the plunger and release the medication.

**Note:** If you accidentally press the auto-release button before you are ready to inject, simply pull the T-bar upward to re-cock the device.

### Step 2:

Press and rotate the tip of the Autoinjector against the skin at the injection site to desensitize the skin and decrease the feeling of needle penetration.

### Step 3:

Press the auto-release button with your thumb (See Figure 6). This will automatically deliver the proper length of the needle into the tissue.

### Step 4:

Reach up with your index finger and completely depress the plunger (See Figure 7). This releases the medication and delivers it into the penis. You can now lift the device away from your skin.

### Step 5:

Immediately apply pressure to the injection site by gently grasping the penis between your thumb and index finger. Maintain firm pressure for at least three minutes.

### Step 6:

Replace the needle cap, pull the T-bar upward and back to remove the syringe, and discard syringe immediately. Make sure that the syringe is not accessible to children once discarded.

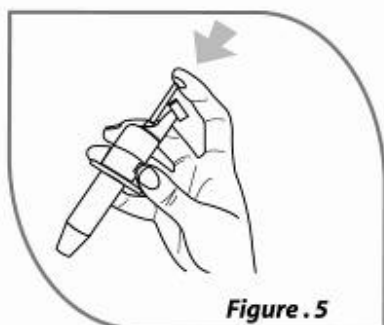


Figure . 5

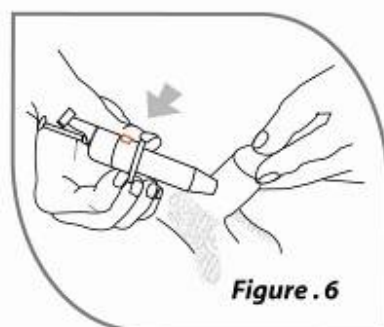


Figure . 6

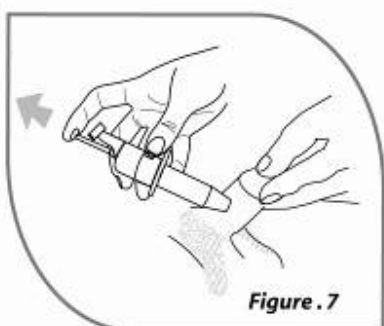


Figure . 7

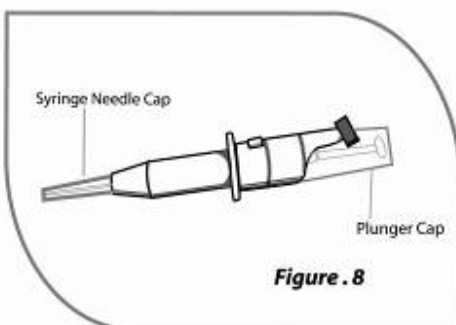


Figure . 8

## CARRYING A PRELOADED SYRINGE

To preload a syringe, fill with prescribed dosage of medication and replace the syringe needle cap. Load the syringe into the Autoinjector and push auto-release button to extend the capped needle. Cover syringe plunger with plunger cap (See Figure 8). Always ensure that your medication is kept cold.

## IMPORTANT REMINDERS

1. Use each syringe once only to reduce the possibility of infection.
2. Always hold the Autoinjector with your thumb resting lightly on the auto-release button. This leaves your index finger free to push in the plunger.
3. Learn to press and rotate the device against the skin just before you press the auto-release button. This action will desensitize the skin and lessen any discomfort from the injection.
4. Always administer penile injections at the **1:00** or **11:00** o'clock position. Rotate your injection sites (Ref. to Figure 4).
5. Consult your doctor immediately if you notice any reactions at the injection site, such as redness, irritation, or swelling or uncontrollable bleeding.
6. Keep syringes away from children before and after use.



**PROLONGED ERECTION!**  
SHOULD YOUR ERECTION LAST LONGER THAN FOUR (4) HOURS YOUR DOCTOR MUST BE CONTACTED IMMEDIATELY ON THE EMERGENCY NUMBER PROVIDED OR CONTACT 083 651 8672 (EMERGENCY ONLY AND IF DOCTOR IS UNAVAILABLE).